

FRONTIER"5513L

High-Speed Micro Centrifuges



Compact Micro Centrifuge with Modern Touch-Wheel Interface for Ideal Routine Centrifugation of Microtubes

Frontier[™] 5513L is a space-saving micro centrifuge capable of spinning 18 x 1.5/2.0 ml microtubes and generating forces up to 15,995 x g for basic molecular biology applications like DNA/RNA purification and extraction. Engineered with convenient operation in mind, Frontier[™] 5513L features a uniquely designed splash-proof front panel with touch-wheel control for quick and easy adjustment.

Unique Features:

- Frontier™ 5513L includes an ergonomically designed 18 x 1.5/2.0 ml microtube rotor that quickly and easily closes with 1/4 of a turn. Adapters allow for spinning of 0.5 ml or 0.2/0.4 ml microtubes.
- Designed for ease of operation, the FC5513L features convenient control on a protective splash-proof interface, which allows for responsive operation even when wearing gloves.
- With a small footprint, Frontier™ 5513L centrifuge is compact enough to save valuable space on the bench, yet large enough to accommodate sufficient number of microtubes ideal for DNA applications.

FRONTIER[™]**55I3L** *High-Speed Micro Centrifuges*

Specifications

Model	FC5513L			
Speed Range	200 rpm – 14 000 rpm; 10 rpm / increments			
Maximum RCF	15 994 x g; 10 x g / increments			
Maximum Capacity (Rotor)	18 x 1.5/2.0 ml			
Temperature Range	Air cooled			
Running Time	Quick Spin, 10 sec to 99 hr 59 min 59 sec, or continuous			
Noise Level	$\leq 57 + 2 dB(A)$			
Acc/Dec Levels	10			
Allowable Density at Maximum Speed	1.2 g/ml			
Allowable Kinetic Energy	1 429 Nm			
AC Power	230 V/50-60 Hz			
Voltage Fluctuation	± 10%			
Current Consumption	0.8 A			
Power Consumption	120 W			
Dimensions (W x D x H)	277 x 356 x 236 mm			
Net Weight (without rotor)	13.5 kg			
Shipping Dimensions (W x D x H)	370 x 470 x 340 mm			
Shipping Weight (without rotor)	15.5 kg			
Environment	For indoor use only			


Rotors & Accessories



Order No.	83041006 Rotor Angle 18 x 1.5/2.0 ml Sealable			
Description				
Compatible with	FC5513			
Max. Speed	14 000 rpm			
Max. RCF	15 994 x g			

Accessories								
Order No.	Description	Vol./Tube	Diameter	Tube Type	No.of Tubes/Rotor			
Without Adapter	Rotor Angle 18 x 1.5 / 2.0 ml Sealable	1.5/2.0 ml	11 mm	Microtube	18			
30130885	Adapter 1 x 0.5 ml D8 mm 6/pk	0.5 ml	8 mm	Microtube	18			
30130884	Adapter 1 x 0.2/0.4 ml D6 mm 6/pk	0.2/0.4 ml	6 mm	Microtube	18			

Outline Dimensions

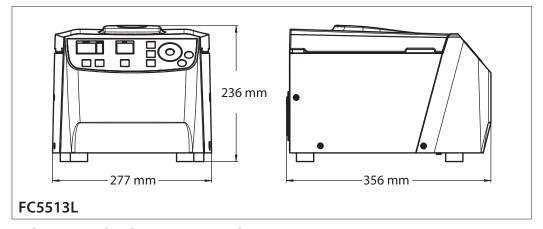

Other Standard Features and Equipment

Ten acceleration and ten deceleration rates, pre-selection of running time from 10 sec to 99 hr 59 min or continuous, audible signal at the end of each run with melody options, quick spin for short runs

Compliance

- Product Safety: IEC/EN 61010-1; IEC/EN 61010-2-020
- Electromagnetic Compatibility: IEC/EN 61326-1 Class B, Basic environments
- Compliance Mark: CE

OHAUS Europe GmbH

e-mail: ssc@ohaus.com Tel: 0041 22 567 53 19 e-mail: tsc@ohaus.com Tel: 0041 22 567 53 20

www.ohaus.com

80776874_A 20230717© Copyright OHAUS Corporation

The management system governing the manufacture of this product is ISO 9001:2015 certified.

